

QC PHANTOM

QUALITY ASSURANCE FOR DR SYSTEMS

The DEL MEDICAL QC Phantom allows comprehensive quality assurance testing for DR detectors and associated software.

The Phantom enables seven quality assurance assessments for monthly and semi-annual evaluation, including exposure linearity and sensitivity, high and low contrast reproducibility, artifact and residual image detection, image resolution, collimator beam alignment, measurement tool accuracy, and display jitter QA testing. The linepair target, 0.10 mm Pb thickness, allows for more comprehensive high contrast resolution assessment.

Key Highlights

Allows Testing of 17" x 17" (43 cm x 43 cm) DR Systems


Includes Test Objects for Seven Quality Assurance Tests

Contains 1 DIN Type 38 Linepair Target

Minimized Lead Content

No Moving Parts

Testing Parameters and QC Report Forms Provided

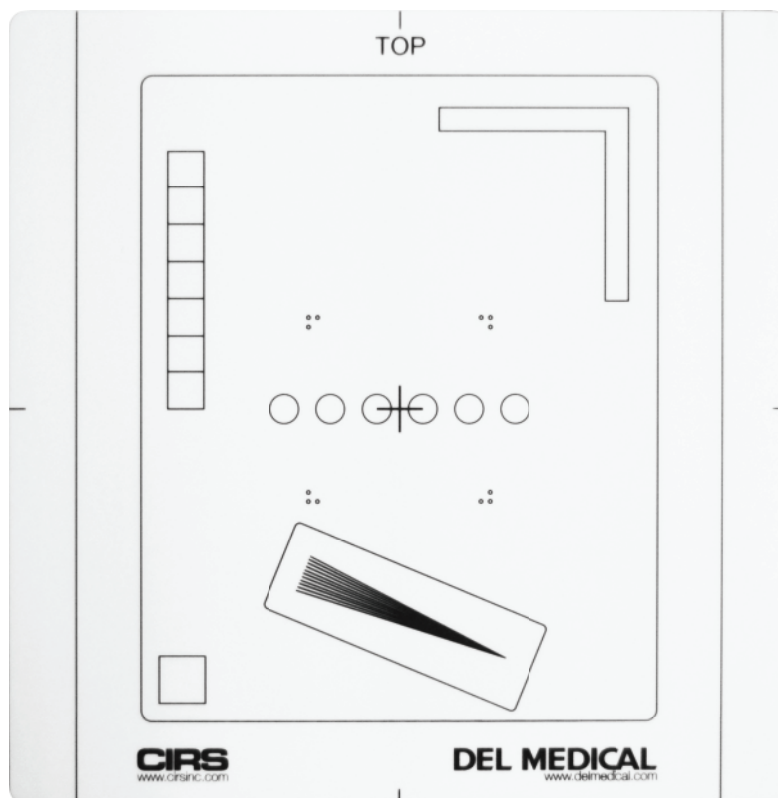

Technical Specifications


* All specifications subject to change without notice

QC Phantom

Dimensions:	17" x 17" x ~ 1/2"
Weight:	8 lbs.
Materials:	PMMA, Copper, Aluminum

Options

Optional accessories:	Padded soft case
-----------------------	------------------

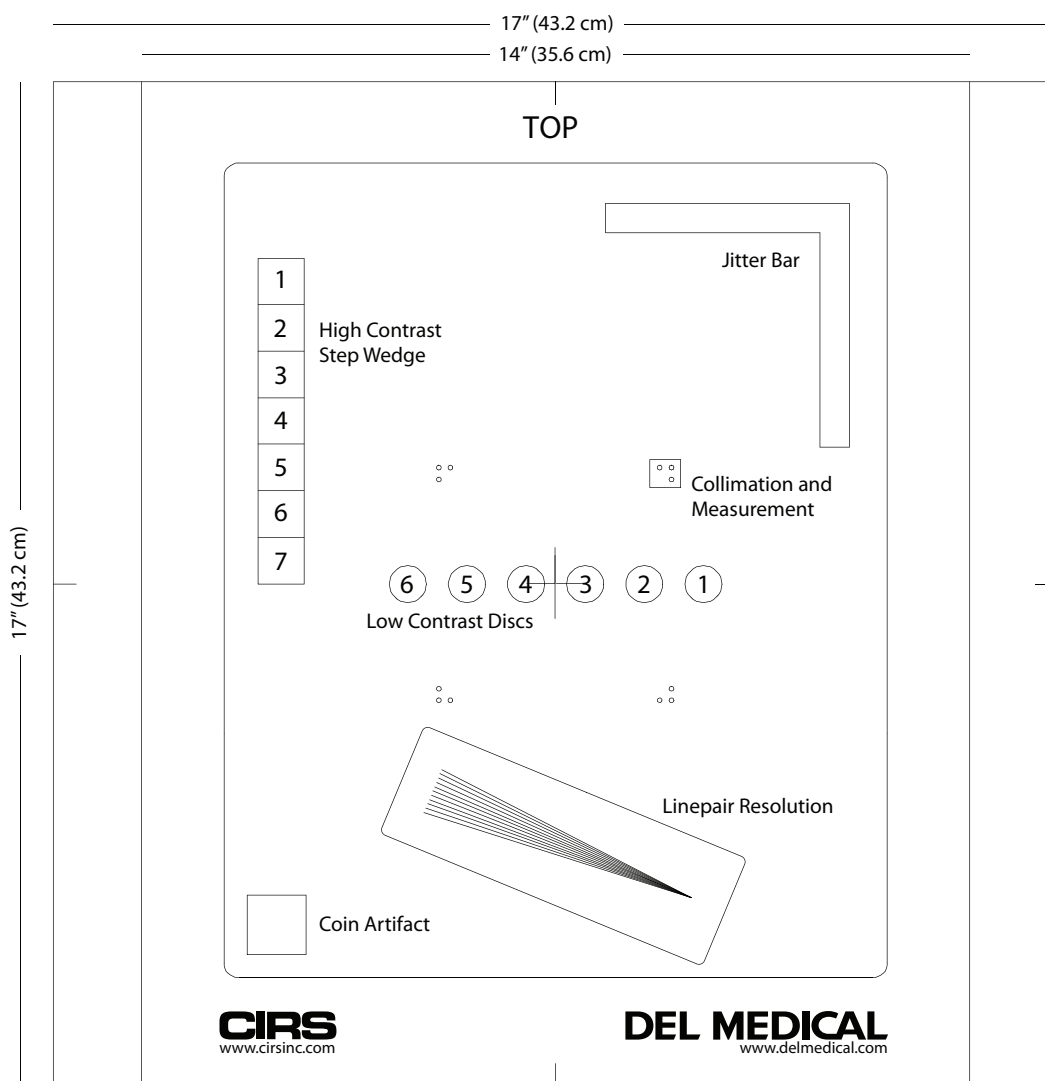

www.delmedical.com

New York
28 Calvert Street
Harrison, NY 10528
(Tel) 800.261.9808
(Fax) 914.835.6111

Illinois
241 Covington Drive
Bloomington, IL 60108
(Tel) 800.800.6006
(Fax) 847.288.7011

DEL MEDICAL

#9000-QCPHANTOM 03.19