

VS50

LONG LENGTH WALL STAND



The DEL MEDICAL VS50 wallstand safely houses the DELWORKS LLI detector for long length and routine upright exams in a fixed location. The DELWORKS LLI allows single exposure long length imaging up to 17 x 42 in (43 x 107 cm) without the need for stitching and is removable for use in other locations.

Key Highlights

Positive Mechanical Lock

Through full range of vertical motion

Full Range of Motion of Active Area of DELWORKS LLI

From 1.5 in (3.8 cm) to 81 in (205.7 cm) above the floor

Fully Concealed Counterbalance System

For easy vertical travel

Wall Hugging Design

Less than 7 in (17.8 cm) from wall to receptor plane

Optional Full-Sized Grid and Convenient Storage Slot When Not Used

Left or Right Load for DELWORKS LLI

Lead Backer for Back Scatter Protection

Acommodates DELWORKS LLI Charging Cable

For continuous charging

No Power Required


Convenient 4-Point Mounting to Wall (No Floor Drilling)

Technical Specifications

* All specifications subject to change without notice

VS50 Wall Stand

Height:	78.36 in (199 cm)
Depth:	7.59 in (19.3 cm)
Receptor Width:	19.12 in (48.56 cm)
Vertical Travel:	37.5 in (95 cm)
Weight (without Grid, Detector, and Lead Backer):	200 lbs (90.7 kg)
Options:	Grid

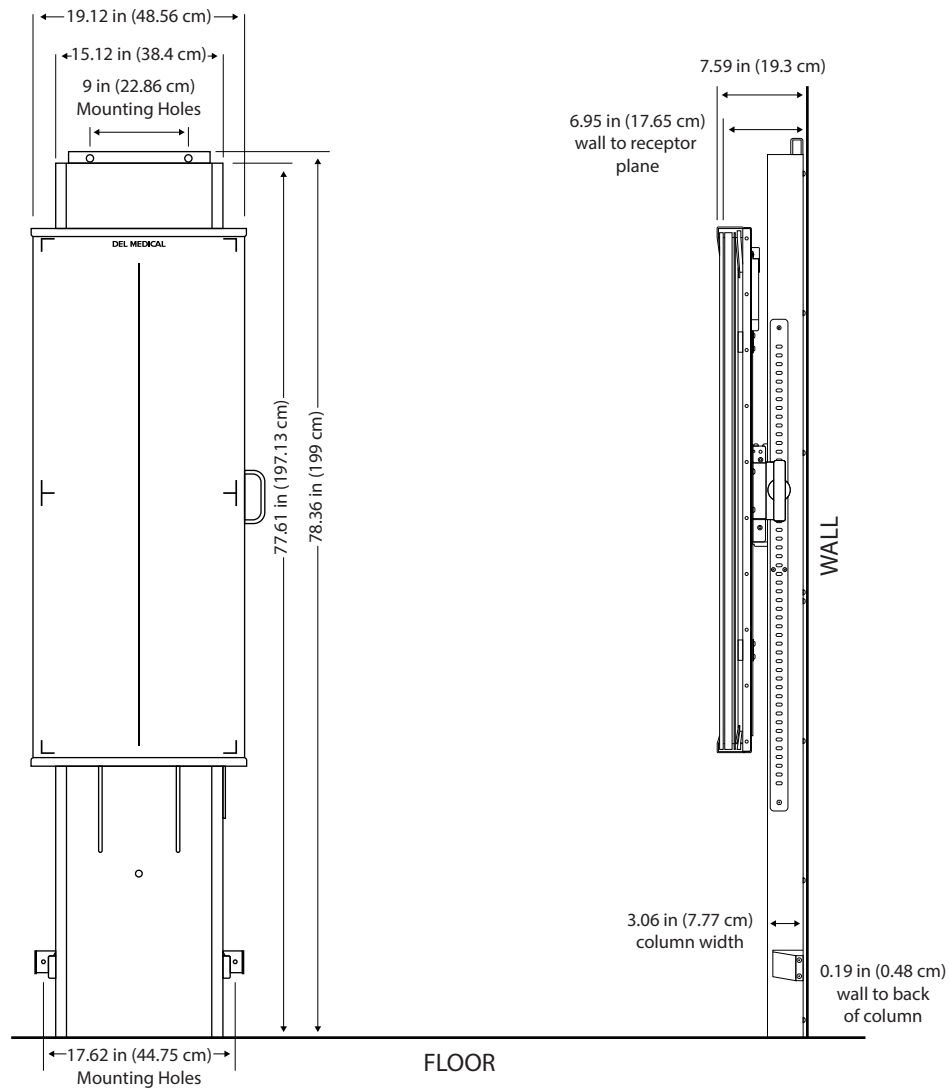

DEL MEDICAL



#9000-VS50-DS Rev. A

New York

28 Calvert Street
 Harrison, NY 10528
 (Tel) 800.261.9808
 (Fax) 914.835.6111

Illinois

241 Covington Drive
 Bloomingdale, IL 60108
 (Tel) 800.800.6006
 (Fax) 847.288.7011

www.delmedical.com