

VS300

PREMIUM WALL STAND

The DEL MEDICAL VS300 wallstand is an easily positioned and secure housing for imaging receptors such as film, CR and all cassette sized digital flat panel detectors. The precise counterbalanced mechanics allows for an effortless and smooth travel of 56.5" from the floor. Soft touch, antimicrobial PA and multi-positional overhead handgrips aide in patient positioning and stability, while the flexible vertical drive positioning handle and conveniently located electromagnetic lock releases reduce technologist wrist strain and fatigue.

Key Highlights

Electromagnetic Locks

- Conveniently placed control buttons are positioned on the side of the receptor cabinet

Ergonomic Release Handle

Side-Mounted and Overhead Patient Hand Grips

- Conveniently placed hand grips to prioritize patient comfort and satisfaction

Vertical Travel Range of 56.5"

Available with 17" x 17" (43 cm x 43 cm) Grid Cabinet or Bucky and a Wide Selection of Grid Options

Cable Concealment System




Technical Specifications

* All specifications subject to change without notice

VS300 Wall Stand

Height (includes Vertical Travel):	83 in (210.8 cm)
Depth:	13.4 in (34 cm)
Depth to Vertical Receptor Plane:	11.6 in (29.5 cm)
Vertical Travel	56.5 in (143.5 cm)
Max Height Center of Receptor	72 in (182.9 cm)
Highest Active Imaging Area from Floor	80.5 in (204.5 cm)
Min Height Center of Receptor	15.5 in (39.4 cm)
Lowest Active Imaging Area from Floor	7 in (17.8 cm)
Cabinet Width	24.6 in (62.5 cm)
Cabinet Height:	26.5 in (67.3 cm)
Weight with Counterweights	200 lb (91 kg)
Weight without Counterweights	N/A
Power Requirement	24 VDC, 1.25 amps max momentary
Front panel aluminum equivalency:	0.4 mm
Options:	Patient overhead handgrip

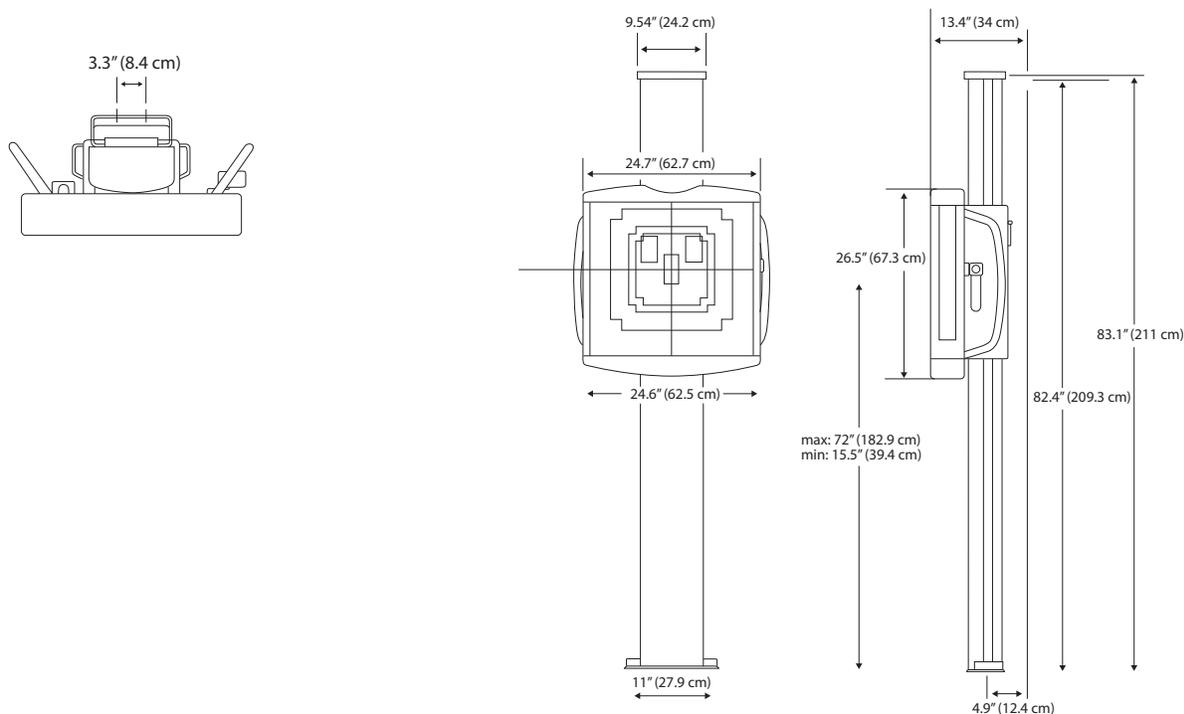

www.delmedical.com

New York
 28 Calvert Street
 Harrison, NY 10528
 (Tel) 800.261.9808
 (Fax) 914.835.6111

Illinois
 241 Covington Drive
 Bloomingdale, IL 60108
 (Tel) 800.800.6006
 (Fax) 847.288.7011

DEL MEDICAL

#9000-VS300 02.19